

fiber<sup>TM</sup>  
D U S T  
THULIUM FIBER LASER



THULIUM LASER

# Fiber DUST™

## THULIUM FIBER LASER




### Dust & Bust

Choose between a fine dusting or fragmentation approach. The use of extreme frequency with very low pulse energy levels broadens treatment options in dusting approach.



### Powerful Popcorning

Powerful and fast reduction of fragment size with popcorning technique.

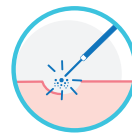

### Soft Tissue Surgery

The output power of 60W allows low power ThuLEP (*thulium laser enucleation of the prostate*).

Lower power emission allows precise ablation and incision of small portions of tissues (such as tumors and strictures).



### Reduced Energy Requirement

Fiber Dust laser is compatible with standard power outlet, with no need of dedicated socket in the OR.

High performances with limited energy consumption!

## Technical Specifications

Wavelength	1.9 $\mu$ m
Emission Modes	Continuous and Pulsed
Average Power	up to 60 W
Frequency	up to 2500 Hz
Energy per Pulse	up to 6J
Cooling	Air cooling
Activation	Double Footswitch

CE  
0123

**Quanta System**  
LASER IN OUR DNA

[www.quantasystem.com](http://www.quantasystem.com)

Quanta System S.p.A - [quanta@quantasystem.com](mailto:quanta@quantasystem.com)

